

1 UNPACK and SET UP

PACKAGE CONTENTS:

- One Karie medication delivery device
- One power adapter: input 110-220V, 50/60Hz, output +5 VDC, 5A
- One set of two security keys
- One Quick Start Guide

MEDICATION CARTRIDGE SUPPLIED SEPARATELY

Read safety information and instructions before use and keep for future reference

SAFETY INFORMATION:


SYMBOL	DESCRIPTION
	<ul style="list-style-type: none">• Caution, sharp edge• Do not touch blade
	<ul style="list-style-type: none">• Do not operate without cartridge

PRODUCT INFORMATION:



Product information label on the bottom of the device


- a** Remove Karie from packaging
Lift straight up out of the foam

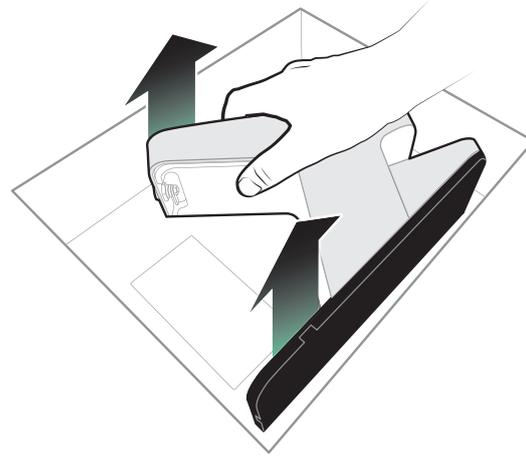


- b** Use Karie upright
Karie is designed to operate upright. DO NOT use Karie on its side, upside down or mounted to a wall



- c** Plug adapter into Karie
Use ONLY the power adapter provided with your device. Fit the power cord into a suitable channel on the bottom of Karie

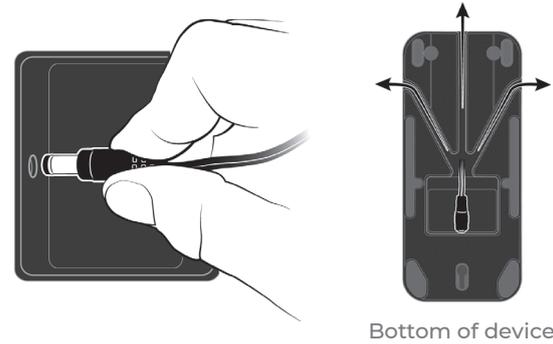

- d** Plug Karie into an electrical outlet
Karie will turn on and connect automatically to a local cellular network. If a 'No Connection' screen appears, see the Karie instruction manual (at kariehealth.com)

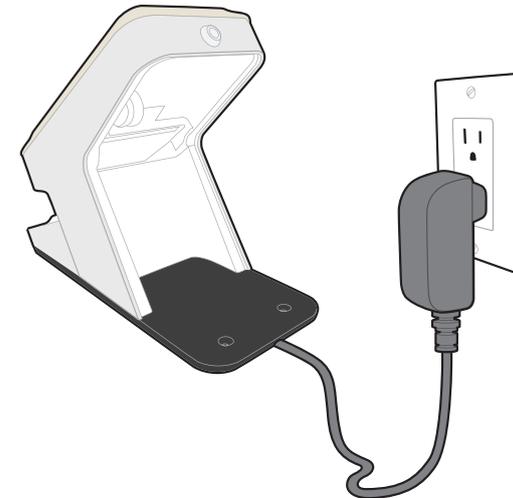

- 2** INSERT CARTRIDGE
Insert the medication cartridge into the back of Karie



Magnets help guide the cartridge into place: you will hear a 'click' when it is fully in place




3 REVIEW TERMS and CONDITIONS and PRIVACY POLICY

To continue using Karie, acknowledge the Terms and Conditions, Privacy Policy and whether to be contacted, by tapping the checkboxes and pressing CONTINUE

After acceptance, Karie will :

- Greet you by name
- Indicate current date & time
- Display your next medication dose time

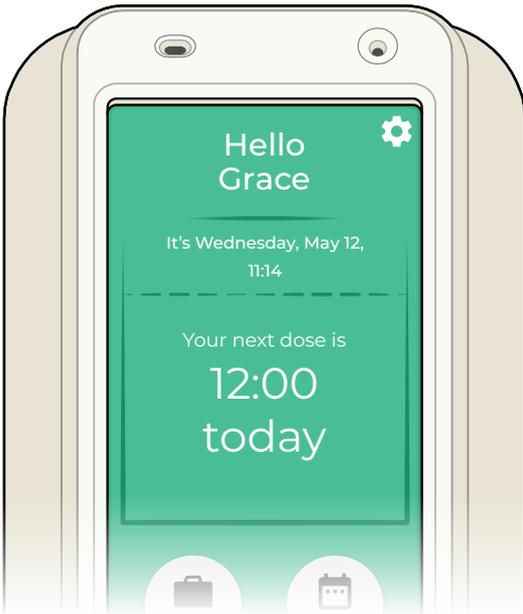


4 CONNECT to Wi-Fi

Karie will connect automatically to a local cellular network. If you live in an area with poor network coverage, follow these steps to connect to Wi-Fi:


a. Press the SETTINGS icon




b. Press CONNECTIVITY

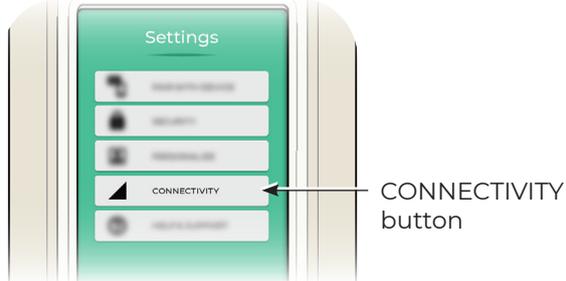

c. Press WI-FI

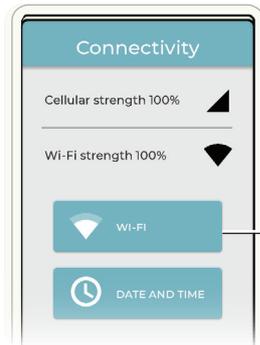

d. Select a network and enter password

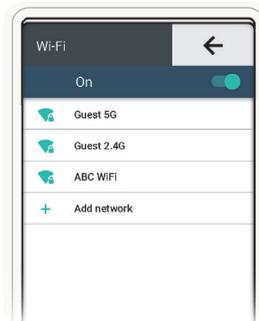

5 TAKE YOUR DOSE


Karie will light up and alert you when it's time to take your dose

Press the Delivery button




Wait for your medication pouch to drop onto the dose tray before taking your medication

Each pouch is pre-cut for easy opening

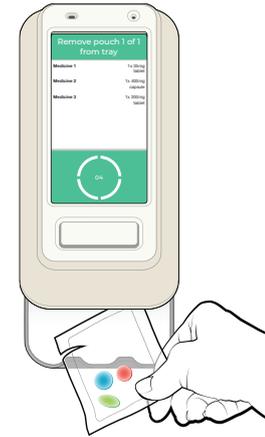

Karie will display your next medication dose time

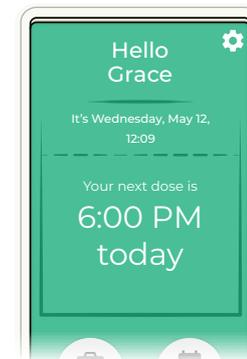

karie™ Quick Start Guide



Visit kariehealth.com/customersupport for the complete Karie Instruction Manual

k karie™ Companion App

- Use the Companion App to pair a caregiver with your Karie device
- You may remove or un-pair a caregiver from your Karie device at any time
- The Karie Companion App is available on the App Store and Google Play



Apple and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries and regions. App Store is a service mark of Apple Inc. Google Play and the Google Play logo are trademarks of Google LLC.